


Suction Machine

SET-UP & RUN



NOTE: Fully charge battery for 17 HOURS. (7314P series only)

1. Securely attach lid to container.
2. Insert container into holder and gently push into place. **NOTE: Do not use excessive force. Pushing container down too hard could cause potential leak or loss of suction.**
3. Attach 4-1/2" (11.43 cm) tubing from filter cartridge to tubing connector on unit.
4. Attach 6 inch suction tubing to the canister lid at the outlet labeled "patient".
5. Ensure power switch is off at this time.
6. Select desired power source. If using AC or DC power, plug the small connector into the DC power input on side of the machine. Plug the other end into an AC wall outlet or DC receptacle.
7. Turn on the suction machine. **NOTE: It's important to make sure you only use grounded outlets for this piece of equipment.**
8. Press and hold the end of the suction tubing with your thumb.
9. Adjust the suction knob to desired level (40-60 for infants, 80-100 for pediatrics, 80-125 for adults).
10. Prior to using, verify that the machine is at the desired suction level and that all connections are secure and without leaks.
11. Attach oral suction catheter/yankauer to end of the suction tubing and suction airway as needed.

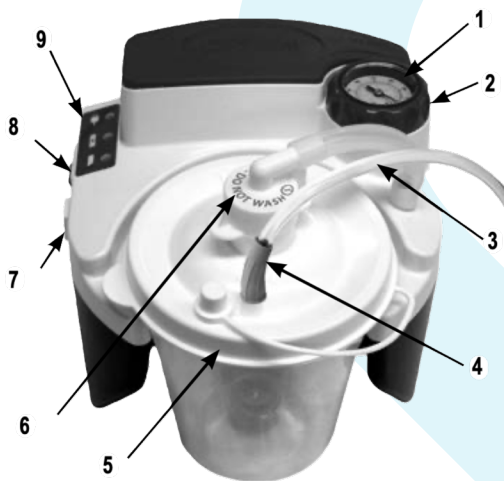

1. Vacuum gauge
2. Vacuum regulator knob
3. 6' patient tubing
4. Patient tubing connector
5. Disposable container with lid (float shut-off of incorporated into lid) and filter cartridge
6. Filter cartridge with 4-1/2" (11.43 cm) tubing
7. DC power input (on side)
8. Power switch
9. LED power lights



Sentido Health

P: 713.955.2123 | F: 281.742.2589

These instructions are not a substitute for the operator's manual.

Credit: DeVilbiss® Healthcare